

# BOTTLE TOP BURETTE



## BOECO DIGITAL CONTINUOUS BOTTLE-TOP BURETTES

- » The digital BOECO Continuous Bottle-top Burette makes manual titration safer, easier, with more reliable results.
- » The DCB burettes can be used for titration agents up to a concentration of 1 Mol/l. During the titration process, the burette fills automatically, eliminating time-consuming refills.
- » Continuous pulse-free operation for faster titrations up to 5 ml per second.
- » Titrant delivery is controlled by turning the wheels, and the amount dispensed is immediately shown on the display - from 10 µl to 999,9 ml, with a velocity of 2,5 ml/turn using model DCB 2500 or 5,0 ml/turn with model DCB 5000.
- » The DCB models eliminate the errors of reading the meniscus and volume graduations of conventional hanging glass burettes.
- » There is no exposed glassware to break and direct connection to the reagent bottle avoids pouring dangerous chemicals.
- » Safety recirculating valve allows for recirculation within the closed system to avoid reagent waste.
- » Angled digital display for easier reading from different heights by different users.
- » Unique telescopic intake tube with length from 210 till 310 mm can be re-adjusted without cutting for different bottle sizes.
- » Dosing spectrum from 10 µl to 999,9 ml.
- » Battery-operated, easily replaceable micro-batteries 1,5 V (LR 03/AAA). One charge lasts a minimum of 500 hours (dispensing 4500 l of reagent).
- » Burette can be rotated 360° on the bottle.
- » Titrating tube removable and adjustable: horizontal 142-220 mm, vertical 10-200 mm.
- » Service friendly modular assembly.
- » Electronic can be completely removed.
- » Low battery indicator.
- » User-serviceable and calibratable to meet ISO and GLP requirements.
- » Optional drying tube protects reagents from moisture or CO<sub>2</sub>.
- » Manufactured under ISO 9001 norms.
- » Adapters are available to fit most bottle neck sizes.

**BOECO Digital Continuous Bottle-Top Burettes, complete with 3 PP adapters 45/38, 38/32, 40/45 (45 thread at the burette), 1 telescopic intake tube (FEP), 1 telescopic titration tube, 2 micro-batteries 1,5 V (LR 03/AAA), without reservoir bottle.**


Code	Model	Accuracy	Coefficient of variation
BOE 9882500	DCB 2500 (2,5 ml/turn)	± 0,2 %	0,1 %
BOE 9885000	DCB 5000 (5,0 ml/turn)	± 0,2 %	0,1 %

### Accessories

Code	Description
BOE 1671095	Drying Tube
BOE 9772050	Amber Glass bottle, 2500 ml
BOE 1671085	Telescopic Intake Tube, 200-350 mm

### Dispenser-thread bottles, with screw cap, amber glass

Code	Volume	Outer diam. (mm)
BOE 9772010	100 ml	28
BOE 9772020	250 ml	32
BOE 9772030	500 ml	32
BOE 9772040	1000 ml	45
BOE 9772050	2500 ml	45

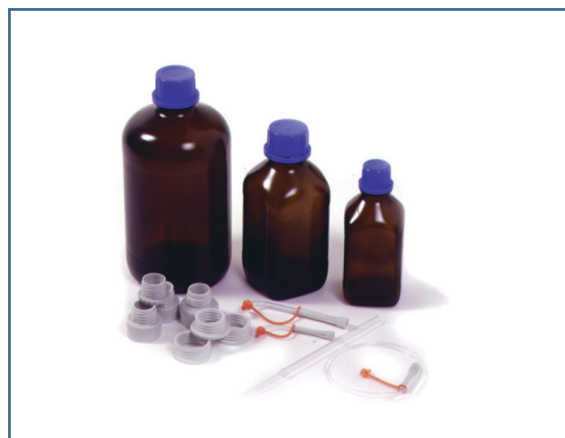

## BOECO

The family-owned company Boeckel + Co. (GmbH + Co.) KG was founded in 1929 and now exists in the fourth generation.

To our customers on all 5 continents we are pleased to introduce our fifth general BOECO catalog with our extended and updated complete line of laboratory equipment and consumables.

We offer high quality at competitive prices by a convenient and cost-effective one-stop supply.

## CONTACT

**Boeckel + Co. (GmbH + Co.) KG**  
**Roedingsmarkt 33**  
**20459 Hamburg**  
**Germany**  
**Phone:** + 49 40 325627-0  
**Fax:** + 49 40 325627-41  
**E-mail:** [info@boeco.com](mailto:info@boeco.com)  
**Internet:** [www.boeco.com](http://www.boeco.com)

**V.A.T. ID-No: DE 118903538**

## BANK

### EURO €

Deutsche Bank AG, Hamburg,  
SWIFT: DEUTDEHHXXX  
IBAN- No. DE 65 2007 0000 0010 5601 00

### USD \$

Deutsche Bank AG, Hamburg,  
SWIFT: DEUTDEHHXXX  
IBAN- No. DE 27 2007 0000 0010 5601 05